


AED Cabinet Specifications



- The measurements provided below are approximate and subject to change without notice.
- Read the "AED Cabinet Installation and Operation Guide" provided with your cabinet prior to installation.
- Check and re-check measurements before cutting holes in walls for recessed cabinets.
- Rough wall opening measurement is the cabinet size. Wall opening will be slightly larger.
- Cabinets equipped with an optional strobe light add approximately 3" to the height.

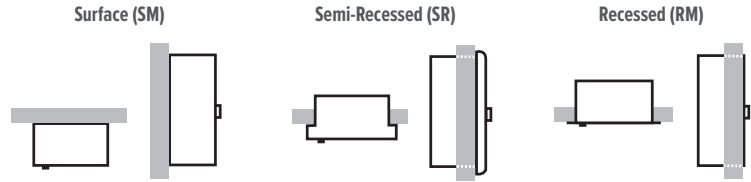

Cabinet Dimensions

Supplier	Style	Mount	Inside Dimensions			Outside Dimensions			Rough Wall Opening			Weight
			Height	Width	Depth	Height	Width	Depth	Height	Width	Depth	
INDOOR CABINETS												
Allied Medical Products	Compact	SM	11.5"	13.875"	6.5"	11.625"	14"	6.75"	-	-	-	8 lbs.
Allied Medical Products	Premium Rounded Corner Compact	SM	11.5"	13.875"	6.5"	11.625"	14"	6.75"	-	-	-	8 lbs.
Allied Medical Products	Standard	SM	14.625"	15.875"	8.25"	14.625"	16"	8.375"	-	-	-	11 lbs.
Allied Medical Products	Premium Rounded Corner Standard	SM	14.625"	15.875"	8.25"	14.625"	16"	8.375"	-	-	-	11 lbs.
Allied Medical Products	Compact	SM	12.875"	13.25"	6.875"	13"	13.375"	7"	-	-	-	10 lbs.
Allied Medical Products	Compact	SR	11.125"	11.75"	3.125"	14.5"	15"	2.875"	11.25"	11.75"	4"	10 lbs.
Allied Medical Products	Compact	RM	11"	11.75"	6.375"	13"	13.5"	.25"	11.25"	11.75"	6.25"	10 lbs.
Allied Medical Products	Standard	SM	17.375"	17.25"	6.875"	17.5"	17.375"	7"	-	-	-	17 lbs.
Allied Medical Products	Standard	SR	14.125"	14.125"	7"	17.5"	17.5"	3.125"	14.25"	14.25"	4"	16 lbs.
Allied Medical Products	Standard	RM	14.125"	14.125"	7"	17.5"	17.5"	1"	14.25"	14.25"	6"	16 lbs.
Allied Medical Products	Tall	SM	30.25"	15.25"	8.25"	30.5"	15.375"	8.5"	-	-	-	23 lbs.
Allied Medical Products	Tall	SR	27"	12"	5.875"	30.5"	15.375"	2.5"	27.125"	12.125"	4"	19 lbs.
Allied Medical Products	Tall	RM	27"	12"	5.875"	30.375"	15.375"	.25"	27.125"	12.125"	6"	18 lbs.
Allied Medical Products	Deep	SM	17.375"	17.375"	8.875"	17.5"	17.5"	9"	-	-	-	18 lbs.
Cardiac Science	Standard	SM	17.375"	17.25"	6.875"	17.5"	17.375"	7"	-	-	-	17 lbs.
Cardiac Science	Standard	SR	14.125"	14.125"	7"	17.5"	17.5"	3.125"	14.25"	14.25"	3.875"	16 lbs.
Cardiac Science	Standard	RM	14.125"	14.125"	7"	17.5"	17.5"	1"	14.25"	14.25"	6"	16 lbs.
Defibtech	Standard	SM	17.875"	17.875"	3.875"	18.25"	18.25"	4"	-	-	-	15 lbs.
HeartStation	Budget	SM	14.875"	13.875"	6.625"	15.375"	14.375"	6.875"	-	-	-	11 lbs.
HeartStation	Compact	SM	12.75"	15.75"	6.25"	16.75"	16"	6.375"	-	-	-	18 lbs.
HeartStation	Compact	SR	13.625"	15.75"	6.125"	19"	18"	2.625"	16.75"	16"	3.625"	21 lbs.
HeartStation	Standard	SM	21.75"	15.75"	5"	22"	16"	6.25"	-	-	-	31 lbs.
HeartStation	Standard	SR	21.75"	15.75"	5"	22"	16"	3"	22"	16"	3"	36 lbs.
Philips Healthcare	Compact	SM	13.375"	16.375"	5.375"	15"	16.5"	6.125"	-	-	-	14 lbs.
Philips Healthcare	Standard	SM	20.875"	15.875"	5.375"	22.25"	16"	6"	-	-	-	27 lbs.
Philips Healthcare	Standard	SR	21.875"	13.875"	5.625"	24.5"	16.5"	2.625"	22"	14.25"	3.75"	30 lbs.
Physio-Control	Standard	SM	17.125"	17.125"	6.125"	17.375"	17.375"	7"	-	-	-	19 lbs.
Physio-Control	Standard	SR	14"	17.125"	6.125"	17.375"	17.375"	3.375"	14.125"	14.125"	3.125"	19 lbs.
Physio-Control	Standard	RM	14"	14"	6.125"	17.375"	17.375"	1.5"	14.125"	14.125"	5.375"	14 lbs.
Physio-Control	Standard	52" Floor	14"	17.25"	7"	-	-	-	-	-	-	32 lbs.
ZOLL Medical	Standard	SM	17.375"	17.25"	6.875" *	17.5"	17.375"	7"	-	-	-	17 lbs.
ZOLL Medical	Standard	SR	14.125"	14.125"	7" *	17.5"	17.5"	3.125"	14.25"	14.25"	3.875"	16 lbs.
ZOLL Medical	Standard	RM	14.125"	14.125"	7" *	17.5"	17.5"	1"	14.25"	14.25"	6"	16 lbs.
ZOLL Medical	Tall	SM	30.25"	15.25"	8.25" *	30.5"	15.375"	8.5"	-	-	-	23 lbs.
ZOLL Medical	Tall	SR	27"	12"	5.875" *	30.5"	15.375"	2.75"	27.125"	12.125"	3.875"	19 lbs.
ZOLL Medical	Tall	RM	27"	12"	5.875" *	30.375"	15.375"	.25"	27.125"	12.125"	6"	18 lbs.
ZOLL Medical	Clear	SM	15.125"	13.375"	7"	16"	16.5"	8.125"	-	-	-	15 lbs.
OUTDOOR CABINETS												
Allied Medical Products	Standard	SM	13.875"	13.875"	9"	16"	16"	9.75"	-	-	-	13 lbs.
WOOD CABINETS												
Allied Medical Products	Compact	SM	12.5"	12.5"	6"	15.625"	16.5"	8.375"	-	-	-	16 lbs.
Allied Medical Products	Standard	SM	15.75"	15.75"	6.5"	18.5"	19.125"	8.625"	-	-	-	24 lbs.

* ZOLL standard and tall metal cabinets include a bracket which reduces depth by up to 1.5"

AED Cabinet Specifications

- Cabinets protect the AED while making it highly visible to potential rescuers.
- Most indoor cabinets are available in powder coated steel. Many are available in stainless steel. Outdoor cabinets are fiberglass.
- AED cabinets are non-locking, allowing quick access to the AED in the event of an emergency.
- Cabinets are available with or without audible alarms and strobe lights. Built-in strobe lights are available for flush mount cabinets equipped with an audible alarm. Satellite strobe lights are available for recessed and semi-recessed mount cabinets equipped with an audible alarm.
- A strobe light is an excellent indicator an alarm is present and a deterrent to tampering. In loud environments, such as manufacturing facilities, it provides a highly visible alert where the audible alarm could go unnoticed.
- A 9V battery provides power for the audible alarm and strobe light. No external A/C power hookup is required. Batteries are supplied with alarms on all cabinets with the exception of Philips brand cabinets.
- Cabinets equipped with an alarm have an “on-off” switch on the cabinet and on the alarm unit.
- Compact cabinets from Allied Medical Products are ideal for these AEDs: Philips Home/OnSite, Philips FR2/FR2+, Philips FRx, Physio-Control CR Plus, and Heartsine PAD.
- Tall style cabinets from Allied Medical Products allow you to add either a small oxygen cylinder or a first aid kit.



Flush-Mount Coated Steel Cabinet with Built-in Strobe



Semi-Recessed Mount Coated Steel Cabinet with Satellite Strobe



Flush-Mount Stainless Steel Cabinet with Built-in Strobe



Flush Mount Fiberglass Cabinet with Built-in Strobe