

ZOLL AED 3[®] Trainer

Training AED and Accessories



Complete ZOLL AED 3 simulation with Real CPR Help real-time CPR feedback

The ZOLL AED 3[®] Trainer is the only automated external defibrillator (AED) training solution that provides Real CPR Help[®], real-time feedback on a trainee's compression quality. The trainer will allow students to learn proper rate and depth of CPR chest compression in a training setting.

The ZOLL AED 3 Trainer is designed to closely mimic the operation of the clinical AED, offering students a more realistic rescue experience. The ZOLL AED 3 Trainer includes:

- An active LCD display with full-colour rescue images that mirror what a student will see on a clinical ZOLL AED 3 device
- Contact sensors in the training pads that detect pad contact with the manikin, eliminating the need for a remote control to support rescue simulations
- An on-device pause button that allows the instructor to pause a training scenario
- Long-lasting consumables, which save time and money
- Easy USB software update/language configuration

2015 AHA and ERC Guidelines Ready

Easy configuration: All configuration and training scenario management is performed on the trainer.

Configurable options include:


- Training scenario
- Operating mode (Semi-automatic or fully automatic)
- Lay-rescuer prompts
- Check-breathing prompts
- Voice-prompt volume
- Language (varies by part number ordered)

The ZOLL AED 3 Trainer includes everything an instructor needs:


- ZOLL AED 3 Trainer device
- ZOLL AED 3 Trainer CPR Uni-padz[®] electrode training harness (minimum of 100 training sessions)
- Trainer electrode gels/liner for ZOLL AED 3 Trainer CPR Uni-padz (minimum of 50 training sessions)
- Four (4) D cell batteries
- Operator's Guide



Real-time CPR feedback with a large 2-colour compression depth indicator

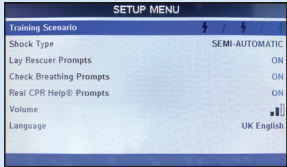

Configure device operation and rescue scenarios directly on the device

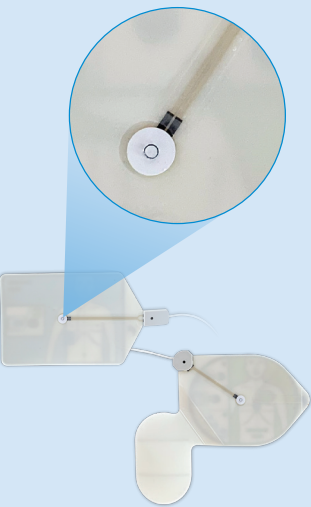

Pad sensors detect contact with the manikin, eliminating the need for a separate remote

ZOLL AED 3 Trainer Specifications

Size: (D x W x H) 5.0 in x 9.3 in x 9.7 in
(12.7 cm x 23.6 cm x 24.7 cm)

Weight: 3 lbs 15oz. (1.8 kg) with batteries

Power: Four (4) D cells (included) for Training AED (Good for more than 16 hours of training)

Rescue Scenarios

Instructors can set up various rescue scenarios directly on the trainer using soft keys.



ZOLL AED 3 Trainer Accessories

Replacements must be ordered separately.



Trainer CPR Uni-padz[®] electrode replacement adhesive gels
Part# 8028-000012 (case of 5)



Trainer CPR Uni-padz electrode replacement liners Part# 8028-000013



Trainer demo manikin
Part# 8000-0835-01



Trainer USB Software Upgrade Kit
Part # 8028-000011-XX



CPR Uni-padz electrode training harness
Part# 8028-000010

To Order Call:

AEDSuperstore[®]

World's Largest Automated External Defibrillator Source

800.544.0048
help@aeds.com

www.AEDSuperstore.com

1800 US Hwy 51 N | Woodruff, WI 54568

Specifications subject to change without notice.

©2019 ZOLL Medical Corporation. All rights reserved. Real CPR Help, Uni-padz, ZOLL, and ZOLL AED 3, are trademarks or registered trademarks of ZOLL Medical Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners. The ZOLL AED 3 defibrillator is not available for sale in the United States.

Printed in U.S.A. MCN PP 1902 0365-05

ZOLL[®]